

# THE BLEEDING STOPS HERE.



QuikClot® Gauze™ Z-Fold Hemostatic Dressing is a 7.5 cm by 3.7 m roll of sterile, white hemostatic gauze. QuikClot® Gauze™ Z-Fold is impregnated with an inert mineral, kaolin, that is both safe and effective in promoting clotting and that does not contain any animal or human proteins. It is packaged in a vacuum packed tearable plastic pouch and is in a z-fold format for easy packing. Also, it contains a x-ray identifiable strip for easy detection.

QuikClot® Gauze™ Z-Fold Hemostatic Dressing is indicated for the temporary control of moderate external bleeding and should be used in conjunction with manual pressure.

# Kaolin is an active agent that creates a strong and stable clot. Our product is:

- Non-invasive
- Not a plug or passive framework
- · Simple, safe and effective
- · Easy to use; Z-Fold format
- X-Ray detectable
- Not known for contraindications
- Vacuum packed

#### **EFFECTIVE**

- Works fast QuikClot® promotes clotting within minutes from application following arterial and venous bleeding even in severe settings
- 5 published independent clinical studies have proven the efficacy and safety of QuikClot® hemostatic agents
- QuikClot Combat Gauze® has been the #1 choice of hemostatic agent for the U.S. Armed Forces since 2008, and continues to be the only hemostatic product approved by the U.S. Military's CoTCCC

#### **EASY TO USE**

- Intuitive simple-to-use dressing format
- · Conforms readily to the wound
- Will not break down or fall apart under pressure

#### SAFE

- No known contraindications and non-invasive
- Inert no human or animal derived proteins no thrombin, fibrinogen or shellfish by-products
- Kaolin, the active clotting agent, has been used safely in cosmetic products for many years
- Independent safety data has been published in the Journal of Trauma (Feb. 2010)

#### **COST EFFECTIVE**


- Safe and less expensive than protein-based products (i.e., thrombin)
- Stops bleeding rapidly, and may reduce the need for more expensive treatments



# THE BLEEDING STOPS HERE.

#### QuikClot® 120 100 No agent added

**Thromboelastography (TEG) Clot Time** 

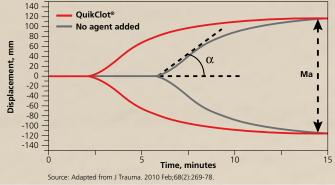

A variety of extensive, independent, pre-clinical testing has been done separately at both the U.S. Army's Institute for Surgical Research (USAISR) and the Naval Medical Research Center (NMRC).

As a result, QuikClot Combat Gauze® hemostatic agent was chosen #1 by the Department of Defense to be carried by all branches of the U.S. Armed Forces, as a first-line treatment for life-threatening hemorrhage.

#### **DIRECTIONS FOR USE**




Remove hemostatic dressing from package



Place the hemostatic dressing into wound and use it to apply manual compression directly over bleeding source. (More than one dressing may be required.)



Continue to apply manual pressure until the bleeding stops.

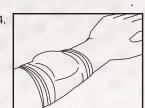

The hemostatic dressing may be left in place for up to 24 hours. To remove, gently remove dressing and thoroughly irrigate wound.



Item #268 QuikClot® Gauze™ Z-Fold **Hemostatic Bandage** Available in -1 unit each

### To place an order:

+1.877.750.0506 Call: Fax: +1.800.343.8656